

Identification of unknown powders

INTRODUCTION

Unknown contaminants such as powders, particles and residues can cause issues in various manufacturing processes. In order to determine the composition of unknown materials, it is often useful to employ several analytical techniques. This paper describes identification of two unknown powders using Raman spectroscopy, combined with two X-ray techniques: X-ray fluorescence spectroscopy (XRF), and X-ray diffraction (XRD). Whereas Raman uses vibrational spectroscopy to determine the bonding environments in organic and inorganic compounds, XRF provides elemental composition. XRD can be utilized to determine phase identification and quantification of crystalline materials based on their unique crystal structure.

ANALYSIS

The two powders of interest in this had a red and black appearance (see Figure 1) and are referred to as Red and Black samples.

Raman spectroscopy is an excellent choice for analyzing unknown powders because it can determine the composition of many inorganic species and is particularly sensitive to iron-based compounds. It is also very applicable to samples that are very black in appearance and may contain large amounts of non-organic carbon. Furthermore, Raman's small spot size makes it

possible to detect rust from particles on the order of a few microns in size.

The Raman spectrum of the black powder is shown in Spectrum 1. The spectrum is consistent with Magnetite (Fe_3O_4) in the spinel structure with a signature peak at 660 cm^{-1} . The red powder Raman results are shown in Spectrum 2. In this case, Hematite (Fe_2O_3) appears to be the main component, although a weak peak at 654 cm^{-1} suggests that some Fe_3O_4 may be present.


To further elucidate the composition, XRF was performed on the two powders. Based on the Raman results, the black powder was likely a magnetite, which means it should have ~ 70 wt% Fe and 30 wt% O. As shown in Table 1, XRF detected 63.5 % Fe and 35% oxygen. The reason for this discrepancy is likely due to the trace contaminants detected in the black powder: Si, Mg, Al, Ti and V were all detected above 0.1 wt%. The combined amount of impurities in the black powder was found to be ~ 1.57 wt%. It should also be noted that whereas Mg, Al and Ti have one oxide form, V and Mn have a range of potential oxide forms with different O content. Thus, the latter are expected to significantly increase the oxygen concentration measured, if they were present.

Similarly, the red powder was found to be primarily hematite with an expected Fe and O content of 72.4 and 27.6 wt%, respectively. The XRF results for the red powder show an even larger increase



Figure 1. Optical images of two unknown powders

Identification of unknown powders

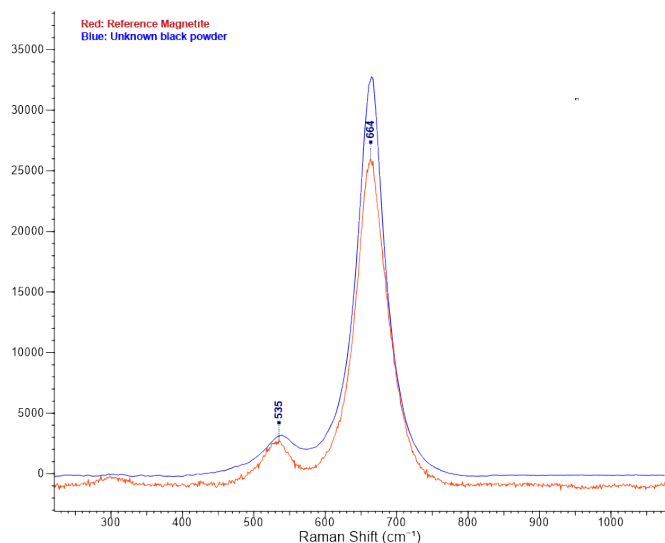


Raman Spectrum 1. Unknown black powder

in oxygen content compared to the expected value. However, this sample only contained 0.68% impurities so the elevated oxygen levels cannot be explained by impurities alone. As suggested by the Raman results, the red powder may not be purely composed of Fe_3O_4 and may contain some Fe_2O_3 , which would increase the oxygen content.

Table 1. XRF results for select major elements

	Black Powder	Fe_2O_3	Red Powder	Fe_3O_4
Fe	63.5	69.9	62.6	72.4
O	35.0	30.1	36.7	27.6
Si	0.44	0	0.24	
Mg	0.27		0.011	
Al	0.21		0.030	
Ti	0.16		0.004	
V	0.12		-	
Ca	0.093			
Na	0.052		0.044	
Mn	0.048		0.15	
Other elements	0.16		0.14	

To determine how much of each iron oxide species is present in the two powders, XRD analysis was performed and the results are shown in Table 2. The reason for discrepancies between the expected and observed Fe:O ratios can be explained by the presence of hematite in the black powder (~7.5 wt%) and a larger amount of magnetite in the red powder (~14 wt%). The detection limit for XRD is roughly ~1 wt%, any contaminants that may be present are likely below this limit.

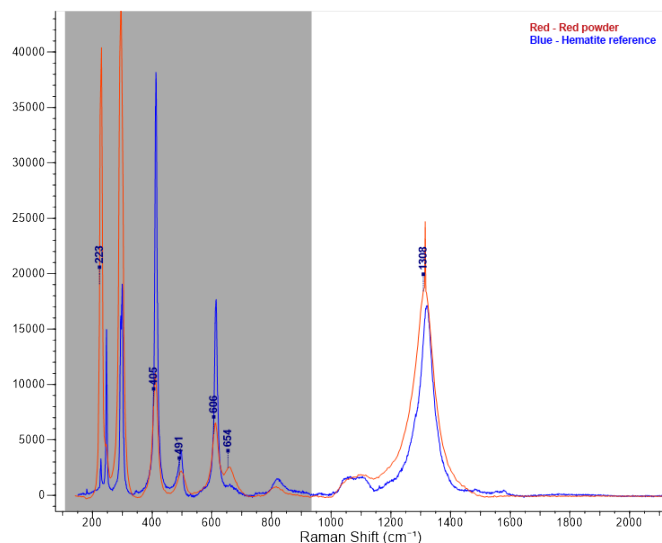

Raman Spectrum 2. Unknown red powder

Table 2. Phase identification and quantitative analysis results

Sample ID	Phases Identified	Concentration (± 5 wt%)
Black Powder	Fe_3O_4 – Magnetite Cubic, S.G.: Fd-3m (227) PDF# 04-009-8424	92.5
	Fe_2O_3 – Hematite Hexagonal, S.G.: R-3c (167) PDF# 01-087-1164	7.5
Red Powder	Fe_2O_3 – Hematite Hexagonal, S.G.: R-3c (167) PDF# 01-087-1164	86.2
	Fe_3O_4 – Magnetite Cubic, S.G.: Fd-3m (227) PDF# 04-009-8424	13.8

Semi-quantitative analysis was performed using WPF (whole pattern fitting), which is a subset of Rietveld Refinement that accounts for all intensity above background. This technique requires that either the structure factors and atomic locations or the reference intensity ratio (a way of comparing the diffracting power of different phases) are known for all phases identified. During this process, the structure factor (related to concentration), lattice parameters (which relate to peak position), peak width and peak shape are refined for each phase to minimize the R value, which is an estimate of the agreement between the model and

Identification of unknown powders

the experimental data over the entire pattern.

The WPF results are shown in Figure 3 and Figure 5 for the Black and Red samples respectively. The concentrations of the phases are listed in Table 1. The R values are ~2% indicating an excellent model fit.

SUMMARY

Thus, while Raman correctly identified the major composition of the two powders, XRF and XRD provided additional details by identifying a mixture of phases in each powder as well as a number of contaminants. This paper demonstrates the synergy of applying multiple analytical techniques to fully understand the chemical makeup of unknown materials.

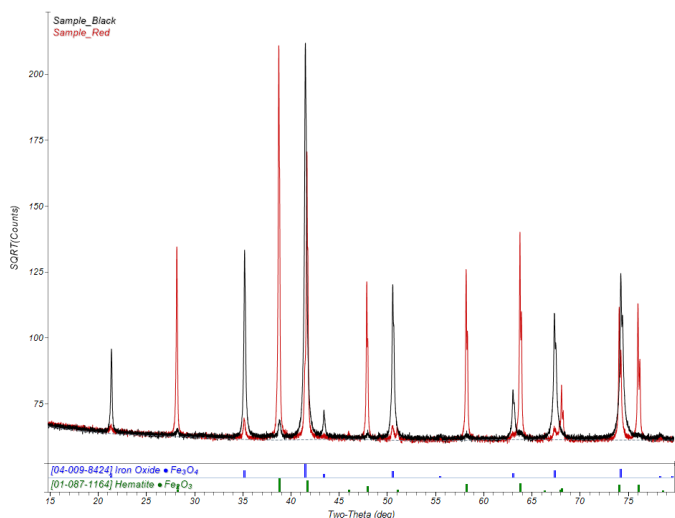


Figure 1: Comparison of diffraction patterns from both powders

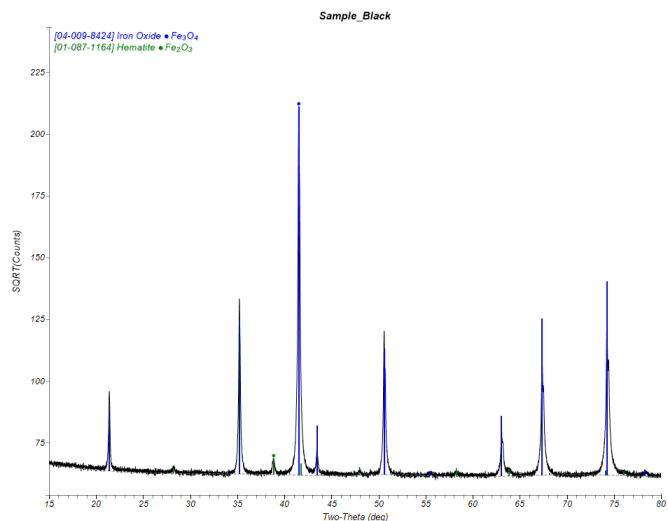


Figure 2: Sample identification for the black powder

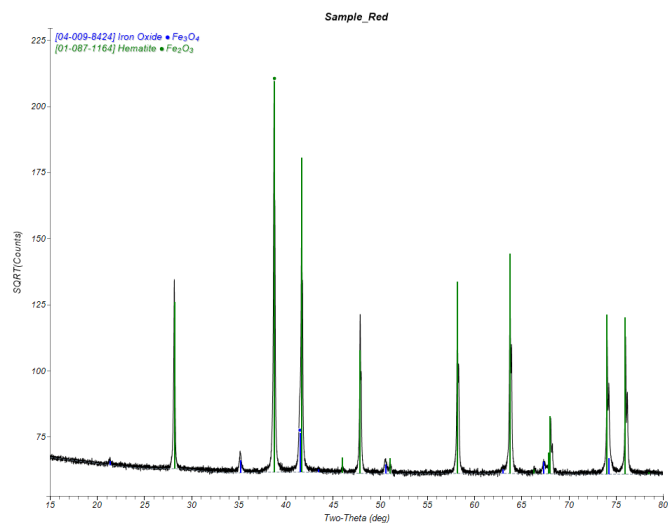


Figure 3: Sample identification for the red powder

Identification of unknown powders

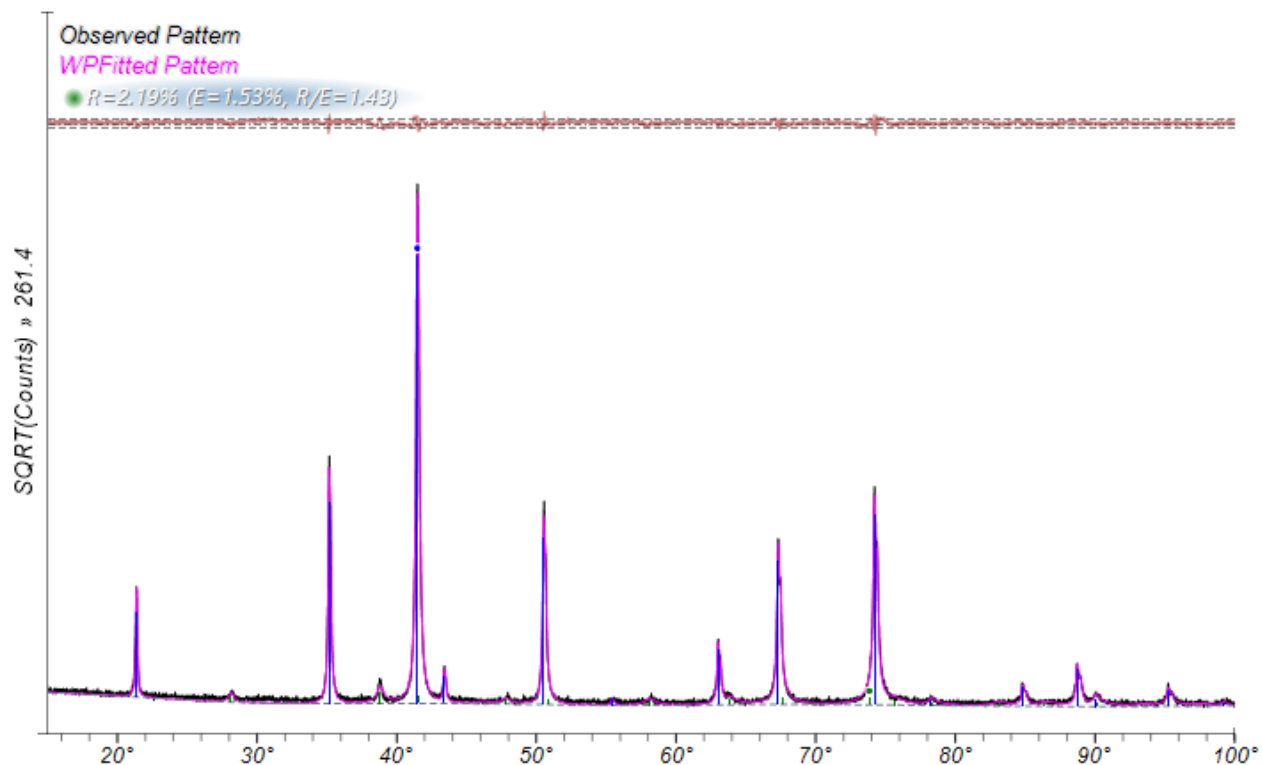
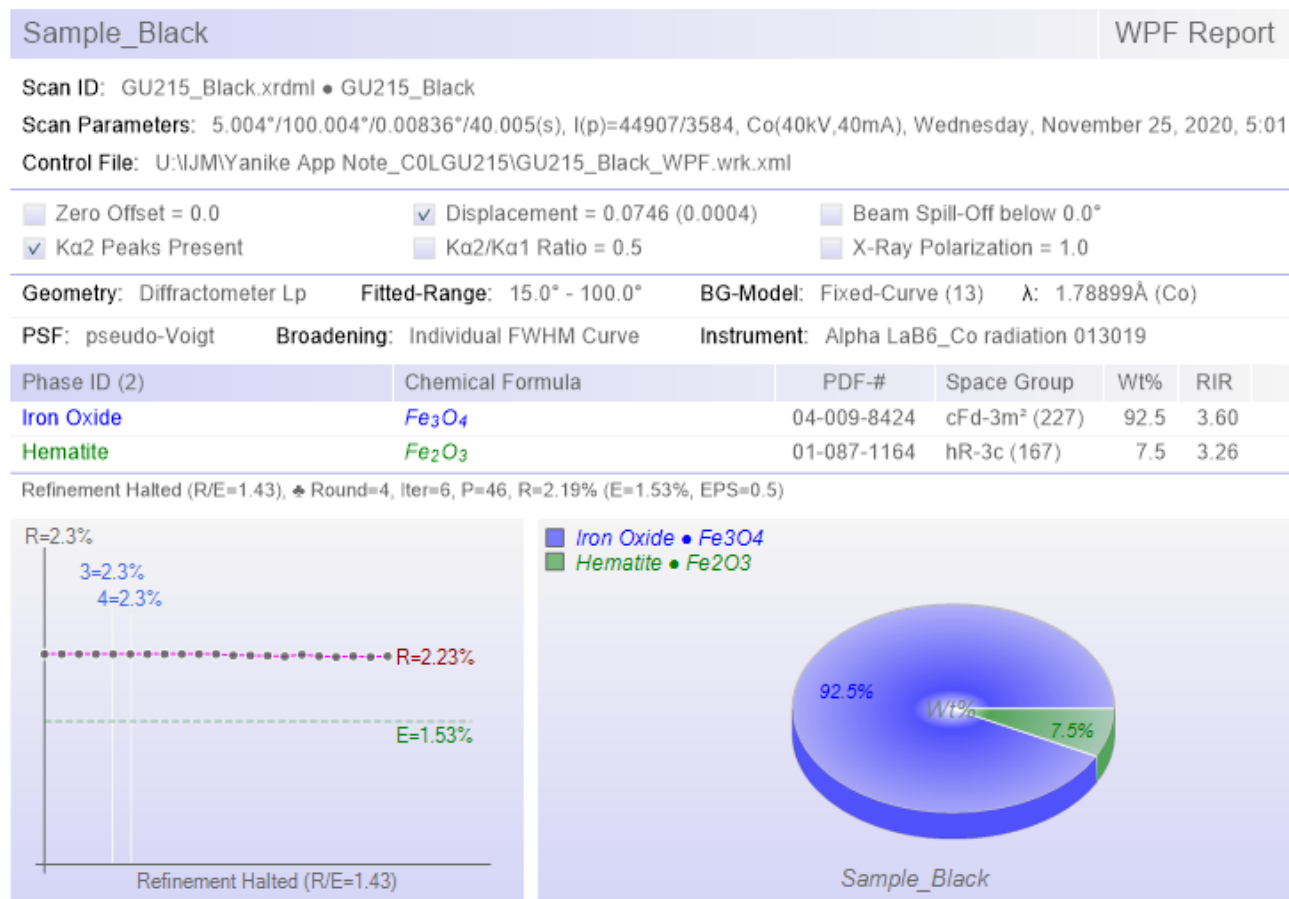


Figure 4: Whole pattern fitting results for black powder

Identification of unknown powders

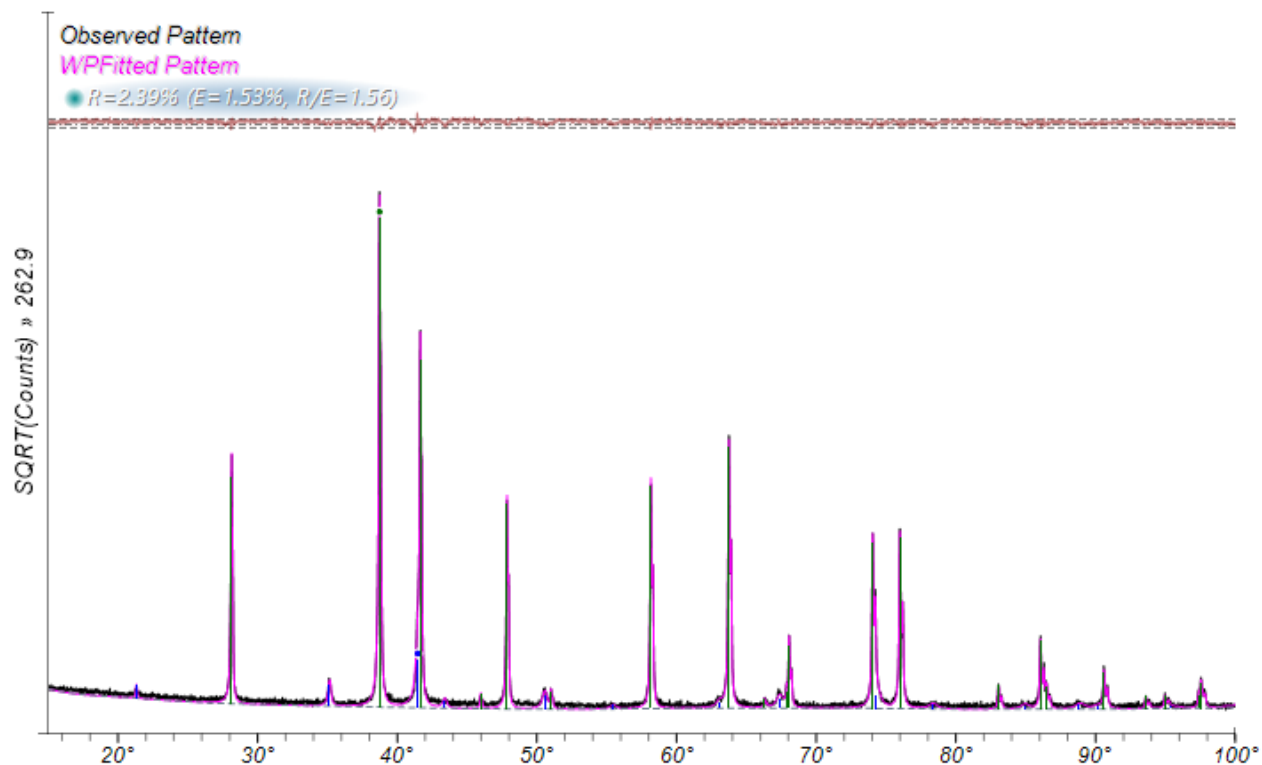
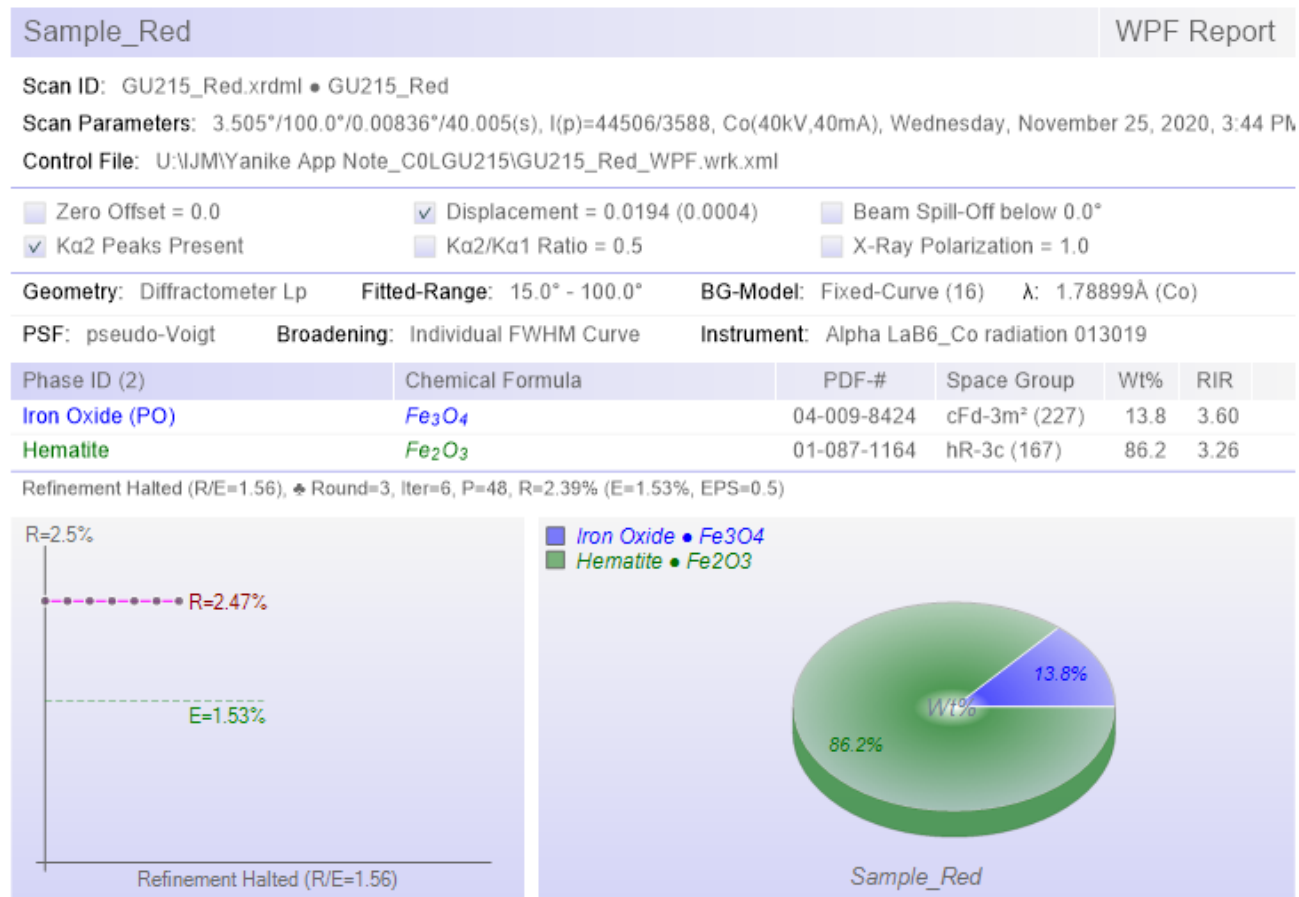


Figure 5: Whole pattern fitting results for red powder